



FAMED

Żywiec

Member of REINSBERG® GROUP



Famed FG-04

Gynaecological chair



The product was made of materials with antimicrobial/antibacterial properties.

www.famed.com.pl

FG-04

Gynaecological chair

The FG-04 gynaecological chair has been designed to support patients during gynaecological and urological procedures. Thanks to standard accessory rails, it can be used for treatments or C-section.

Chair can be fitted with different accessories, many of them are presented in the accessories for operating tables catalogue.

Ease of adjustments

- For the ease of use you can use electric control of: height adjustment, independent seat inclination and back segment inclination. It is also possible to achieve both Trendelenburg and reverse-Trendelenburg.
- The back segment inclination is regulated simultaneously with seat segment inclination- a constant back pressure prevention.
- Chair can be regulated from remote control or foot controller.

Structure of the chair

- The chair is mobile - you can easily change its location in the office. It can be locked in place using one single lever that do not disturb the doctor (it is located at the back of the chair).
- Profiled front of the base makes it easy to climb into the chair, it also allows the doctor to rest his feet against the base.
- All moving parts and main construction of the chair is covered by ABS - it is manufactured with pSilver™ antibacterial additives.
- Stainless steel accessory rails are available on both seat section and backrest.
- Integrated rails for stainless steel bowl to catch fluids.

Easy to clean, high quality upholstery, different colors available

Different leg holders or arm supports available

Stainless steel bowl, there is also a bowl that can be mounted on te rails

Moving elements and actuators fully covered by antibacterial ABS, very easy to clean

Profiled base - easy entry, possibility to rest doctor's feet on it

Upholstery

- The chair is equipped with seamless upholstery that is antibacterial (pSilver™), antiallergic, phthalate-free. It is produced in Poland.
- To fit the modern medical environment, you can choose variety of colours of the upholstery.

Advantages of the FG-04

- Wide, comfortable seat and easy access.
- Easy positioning due to 3 electric actuators.
- Wide range of upholstery colours.
- Seamless upholstery and ABS covers facilitate keeping the chair clean.
- Applied materials are resistant to disinfecting agents.
- Wide range of equipment allows the user to choose optimal version of the chair.
- Option to use it as a treatment table.
- Wheels facilitate displacing the chair in the cabinet.
- High SWL enables you to treat obese patients.
- Memory functions in the remote.



Integrated paper roll holder, fully hidden

Contoured back section to allow easy access for anesthesiologist



Advanced wired remote control with memory functions

Accessory rails compatible with OR-specific accessories

Foot controls



Wide seat for comfort

Positioning from remote

High quality upholstery

Mobile base with lock

High quality actuators

OR accessories can be used

Fully hidden couch cover

Easy disinfection

FG-04 with hand rests and knee rests (WF-42.0)



Standard configuration


- Foot controller.
- Wired remote control with 2 memory functions and custom memory button.
- Paper roll holder, completely hidden to ease the disinfection process.
- One piece of couch cover paper roll.
- 4 stainless steel rails that can fit surgery-specific accessories.

Additional accessories and options

- Bowls for fluids [WF-52.2] - 4 liters in capacity.
- Foot stirrups [WF-23.0] - mounted on leg rest (no way to use them without leg rest, which is a part of WF-42 set).
- Leg rests and arm supports [WF-42] - leg rests are mounted to the accessory rail using clamp, arm supports are mounted to the seat section.
- Foot rests for lithotomy positioning, with abduction possibility. PU covering of the stainless steel rods next to patient body, antibacterial footprints [WF-36.0].
- Bowl for fluids [WF-52.2].
- Foot stool with upholstery [WF-03.3].
- IV pole with WS-17.6 clamp (with angle adjustment) [WK-01.5 W2].
- IV pole with WS-16.5 clamp (without angle adjustment) [WK-01.5 W1].
- Chair can be equipped with battery [version no. FG-04 MW8] or without [version no. FG-04 MW5].



Wired remote control

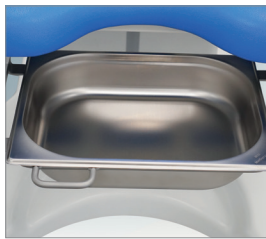

Urological bowl WF-52.2 [option]



Stainless steel foot stirrups WF-23.0 - mounted to leg rest [option]



Leg rest (pair, mounted on rail) and arm support (pair) (WF-42) [option]



Foot rests with abduction and antibacterial footprints (WF-36.0)



Foot stool with upholstery WF-03.3 [option]



IV pole, WK-01.5 W1



Side rails for accessories in the seat and backrest.



Central locking system



Paper roll storage



FG-04 with leg plate extension (WF-22.3), set in a Trendelenburg tilt.


Chair can act as a delivery table, C-section table and simple examination and treatment unit.

Accessory rails can adapt OR-specific accessories to support the patient during treatment.



Upholstery

Antibacterial upholstery, CRIB 5

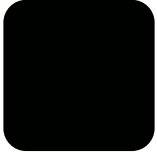

Black
H1/Ab1



Heather
H4/Ab4



Aquamarine
H6/Ab6



Tangerine
H8/Ab8

Antibacterial upholstery, phtalate-free



Yellow
1F



Light blue
2F



Blue
3F



Turquoise
4F



Light grey
5F



Grey
6F



Black
7F




Red
8F


Electric functions and ROM




Functions executed by actuators and their range of motion

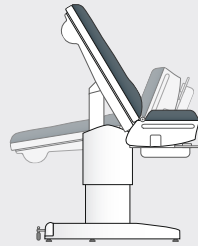


-10° to +65°

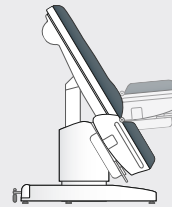


-65° to +10°

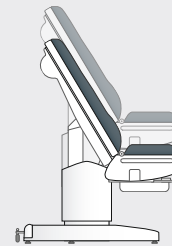

560 to 860 mm
with seat in
horizontal position

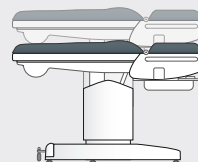

780 to 1080 mm
with seat and
back section in
horizontal position

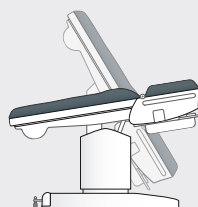

65°

Technical data

Total seat and back section width	745 mm
Width (seat section / back section)	570 mm / 565 mm
Height adjustment in chair position (in horizontal position)	560 - 860 mm (780 - 1080 mm)
Back rest inclination	-10° to +65°
Seat segment inclination	-65° to +10°
Trendelenburg / reverse-Trendelenburg position	10° / 65°
Power supply	230/127/110/100 V ~, 50/60 Hz
Maximum power consumption	300 VA
Electric protection class / application part	II / B
Environment influence protection degree	IP X4
SWL	200 kg

Edition 05/2024/02/4. Famed Żywiec reserves the right to modify the product and specification as part of technical progress. All illustrations and photos used in this material are for illustrative purposes only and may not reflect the finished product. The people shown in the photos are not medical professionals. They are models. The device presented in the catalog is intended for use in health care facilities by authorized persons after having read the manual (IFU).

